

X-Ray Illuminators

HXR400

Healthcare Lighting's HXR400 X-ray Illuminators are mounted flush with the wall and provide a clean appearance to the room. The HXR400 with (4) lamps per section is standard. Choose from our extensive options and build an illuminator that meets your individual needs.

(FLUSH MOUNT ONLY)




FEATURES

Specifications

(Specifications subject to change, consult web site for any revisions)

- Construction:** Continuous housing and individual reflectors are 20GA corrosion resistant, electro-galvanized steel; all corners are welded for added strength; units longer than (6) horizontal sections are separated into (2) pieces for shipment.
- Lenses:** Continuous 1/8 inch white translucent Plexiglas; internal section dividers are removable for panoramic film viewing; units longer than (6) horizontal sections are separated into (2) pieces for shipment.
- Electrical:** Flicker free electronic ballasts are standard; double insulated low leakage ballast is available for areas with EMI sensitivity; optional electronic dimming to detect subtle differences in tissue and bone densities.
- Switching:**
- One durable, oversized rocker switch is standard per section.
 - Film activated switches offer instant-on when film is inserted.
 - Console switching conveniently places all switches at the bottom center of the housing with an additional, master-switch that controls the entire unit (recommended for units greater than (4) sections).
 - Desk-mounted Remote Switches are available for greater operating ease. Individual section switches are included as well as master switches that control the entire fixture.
- Lamps:** (4) 15 watt T8 daylight fluorescent lamps per section are standard. (See below for optional lamp configurations; lamps included).
- Film Retention:** Exclusive adjustable roller film grips.
- Mounting:** Available in flush mount only; *Not available in surface mount.*
- Standard Finish:** Trim face is powder coated white.
- Listings:** UL Listed.

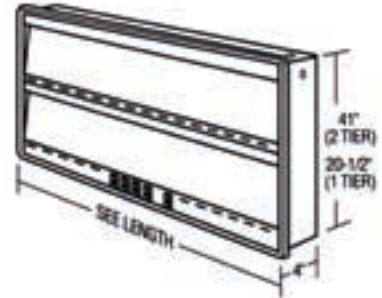

HXR404-S-2-CSW shown above

Note: For Flush Mounted Units

Add 3/8" overall to fixture dimensions for flush fixture rough-in openings. Flush models have 3/4" flange to cover rough-in gap.



ORDERING INFORMATION

Series

Flush Mount ⁽¹⁾	Nom. Hsg. Dim. (L x W x H)	Lamps
HXR401-R-1 (single tier).....	1-Section (14-5/8" x 4" x 20-1/2").....	(4) 15w T8 (lamps included)
HXR402-R-1 (single tier).....	2-Section (28-3/8" x 4" x 20-1/2").....	(8) 15w T8 (lamps included)
HXR403-R-1 (single tier).....	3-Section (42-7/8" x 4" x 20-1/2").....	(12) 15w T8 (lamps included)
HXR404-R-1 (single tier).....	4-Section (57" x 4" x 20-1/2").....	(16) 15w T8 (lamps included)
HXR405-R-1 (single tier).....	5-Section (71-1/8" x 4" x 20-1/2").....	(20) 15w T8 (lamps included)
HXR406-R-1 (single tier).....	6-Section (85-1/4" x 4" x 20-1/2").....	(24) 15w T8 (lamps included)
HXR407-R-1 (single tier).....	7-Section (99-3/8" x 4" x 20-1/2").....	(28) 15w T8 (lamps included)
HXR408-R-1 (single tier).....	8-Section (113-1/2" x 4" x 20-1/2").....	(32) 15w T8 (lamps included)
HXR401-R-2 (double tier).....	2-Section (14-5/8" x 4" x 41").....	(8) 15w T8 (lamps included)
HXR402-R-2 (double tier).....	4-Section (28-3/8" x 4" x 41").....	(16) 15w T8 (lamps included)
HXR403-R-2 (double tier).....	6-Section (42-7/8" x 4" x 41").....	(24) 15w T8 (lamps included)
HXR404-R-2 (double tier).....	8-Section (57" x 4" x 41").....	(32) 15w T8 (lamps included)

Options

- MSWFilm Activated Switching (per section)
- 2CSW2 Section console switching
- 3CSW3 Section console switching
- 4CSW4 Section console switching
- 5CSW5 Section console switching
- 6CSW6 Section console switching
- 7CSW7 Section console switching
- 8CSW8 Section console switching
- 2/2CSW2/2 Section console switching
- 4/4CSW4/4 Section console switching
- 6/6CSW6/6 Section console switching
- 8/8CSW8/8 Section console switching
- RTSRemote Switching (per section)
- ELDBi-Level electronic dimming ballast (per section)
- LLBDouble insulated, low leakage ballast (per section)

Ordering Notes:

- (1) The HXR400 Series is available in flush mount only. Not available in surface mount.

