

# Major O.R.

The Visionary® supplies light necessary for any major surgical procedure and reduces the risk of errors. Each head provides 8,300 foot candles at one meter. With a double head mount, Visionary® incorporates (8) lamps total on (2) heads with no shadows.

Offering outstanding illumination and an array of superior features, the Visionary® is ahead in innovation and value. It enables a facility to upgrade to a surgical lighting system that will satisfy even the most demanding surgeons.

*(Formerly known as Odyssey)*



## FEATURES


### Single Ceiling Mount



### Double Ceiling Mount



### Wall Dimmer Switch

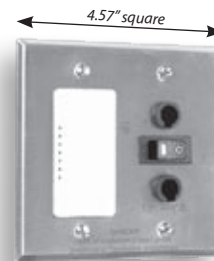

## Specifications

*(Specifications subject to change, consult web site for any revisions)*

- Housing: Injection molded polycarbonate exam housing; yoke and arms are welded steel.
- Rotation: 360° Limitless mounting system rotation around vertical axis *(single ceiling models only)*.
- Switching: **Wall mounted dimmer switch is included; 0-100% continuous dimming.**
- Lamp: Single Visionary® head includes (4) 35W, BA15d base; 24 volts, CC-8 halogen lamps wired in parallel for continuous operation; rated for 1000 hours; 4200°K, 98 CRI; output 8,300 foot candles at 1 meter *(lamps included)*.
- Focusing: Adjustable focus with rotation of center autoclavable handle.  
Light field diameter: 8 inches to 10 inches.  
Depth of illumination: 22 inches.
- Mounting: Available in single and double ceiling systems.  
Ceiling mounted models height restrictions: 8.5 feet to 10 feet.
- Reach & Sweep: Ceiling mount horizontal radius: 54 inches.  
Ceiling mount vertical sweep: 47 inches.
- Finish: Housing, yoke, arms and handle are baked enamel powder coated light grey.
- Handle: One autoclavable handle comes standard with every unit; disposable handle covers are available.
- CA OSHPD: California Office of Statewide Health Planning and Development: Mounting system meets California seismic codes.
- Manufacturer: Visionary® is manufactured by Burton Medical, Inc. Healthcare Lighting® is a distributor for this product.
- Listings: UL Listed: 60601-1.

*Lamp maintenance is a snap. Just push and twist the rotating lens cover, remove the old lamp and insert the replacement into the socket.*



## ORDERING INFORMATION

*Example: HOR301-115-9014*

Series <sup>(1)</sup>

..... Description .....

- HOR300-115** Single Visionary® head (Ceiling Mount-120v) (4) 35w, 24v quartz halogen lamps (included)
- HOR301-115** Double Visionary® head (Ceiling Mount-120v) (8) 35w, 24v quartz halogen lamps (included)

Options

- 8100** Box of 25 SteriHandle™ disposable covers
- 9014** Replacement autoclavable handle
- 2000** Replacement 35w, 24v quartz halogen lamps (pack of 3)

**Ordering Notes:** (1) Product available in 230v (International Sales Only)

