





SOLUME XT MOTORIZED EXAM

Dimensional Data

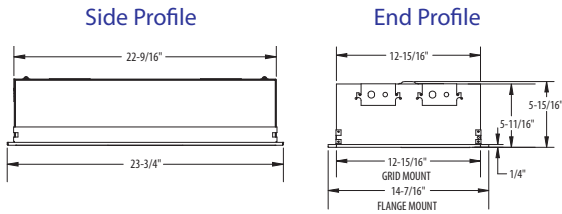

Range of Motion

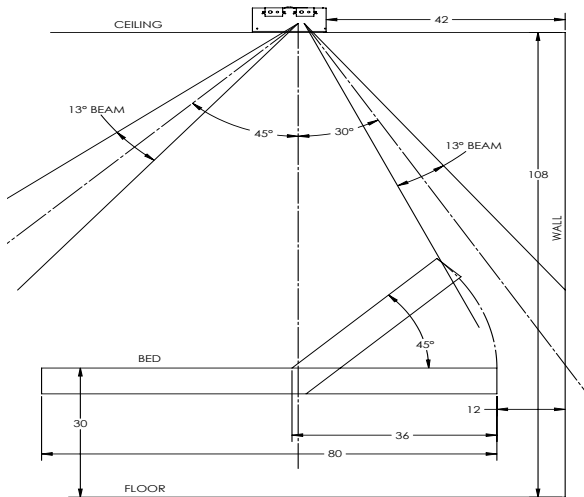

Specifications subject to change. Please visit our web site for updated information.

Specifications

- Description:** Recessed, motorized exam light.
- Range of Motion:** Exam head rotates from -30° to 45°.
- Applications:** General patient room and Med/Surg patient room.
- Construction:** 18GA welded housing and exam head.
- Electrical:** 120 or 277 volt.
- Lamp Shield:** Examination head incorporates a tempered glass lens to shield lamps.

Exam:

MA:

Motorized Examination System includes:

- Healthcare Lighting's Class 2, four Circuit Low Voltage Controller which provides interface to Class 2 Wall Switch for control of motorized examination head. Class 2 Wall Switch fits a standard 2 gang electrical box (by others).
- (3) 75w CRI 100, MR16 lamps (included) with heat dissipating dichroic filters. 13° beam, 4,000 hour lamp life.

Chart Light:

Optional (1) 3w Cool White LED or (1) 10w MR11.

Finishes:

Polyester powder coated after fabrication with either flat white or designer white antimicrobial finish.

Listings:

UL and C-UL listed.

Warranty:

5-year limited warranty. Made in the USA by Healthcare Lighting's ISO certified facility in Fairview, PA.

SOLUME XT



Ordering Information

Example: HMS601-MA-F-120-LEDCR-AM

1	2	3	4	5	6
HMS601					
SERIES	EXAM TYPE	MOUNTING	VOLTAGE	OPTIONS	FINISHES
HMS601	MA Motorized	G T-Grid F Flange	120 277	LEDCR LED Chart Light on right LEDCL LED Chart Light on left CR11 MR11 Chart Light on right CL11 MR11 Chart Light on left EC Remote pendant exam control cord with holder bracket (MA only)	FW Flat white finish AM Designer white antimicrobial finish